

We recommend the use of our duplicating silicones, 1:1.
Available as:

- 1 kg, component A / B
- 6 kg, component A / B
- 25 kg, component A / B



WaSil soft,
hardness Shore A: 16-18,
6 kg cask



WaSil speed,
hardness Shore A: > 24,
6 kg cask




Sidomix with 6 kg cask



Single-handed operation



Material container, 2 x 2 L,
incl. in delivery

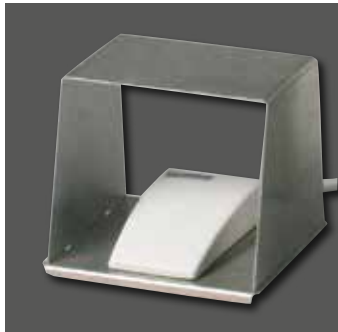

Foot-operated switch,
available as option

Wassermann is a medium-sized, owner-operated company with head office in Hamburg. We design and produce high-quality technical devices in Germany. Our main users are dental laboratories, dental clinics and various industrial areas worldwide.

Why a Wassermann product?

Wassermann products are manufactured from high-quality materials and the proportion of plastics used is intentionally kept as low as possible. Our devices are known for their suitability for daily use, functionality and longevity. As a result, the users place high levels of trust in the Wassermann brand. We feel it is our duty to support our customers on a daily basis when it comes to safety, service and productivity. Technological perfection, optimum quality and especially sustainability are the benchmarks we set. Thanks to the combination of these values, the family-owned company, since its establishment in 1927, has gained worldwide recognition as a manufacturer of high-quality products.



Further information are available on our website at
www.wassermann.hamburg



/WassermannDentalmaschinen

Your Dental Dealer _____

Sidomix



Dosing and Mixing Unit



Wassermann Dental-Maschinen GmbH
Rudorffweg 15-17 · 21031 Hamburg · Germany
Phone: +49 (0)40 730 926-20/ -24 · Fax: +49 (0)40 730 37 24
export@wassermann-dental.com · www.wassermann.hamburg
www.facebook.com/WassermannDentalmaschinen

Dosing and Mixing Unit

Sidomix

The Sidomix is a dosing and mixing unit for addition-vulcanizing duplicating silicones of low viscosity. A regulable bypass circuit ensures an optimum mixing ratio. The mixing ratio can be adjusted very easily and quickly by the user with the help of a convenient quantity flow control. The desired components are mixed 1:1 and dosed correctly – automatically, totally homogenous and fully bubble-free.



Convenient: The containers for the A and B components can be fitted very easy to the unit by the quick couplings

The mixing head with separate outlets makes mixing of components impossible – even when attaching and removing the mixing cannulas

Material and time saving with loss-free filling of several flasks with one mixing cannula

High quality construction, plastic powder coated housing

Major Features

- **Included in delivery:**
Set material cask with quick fitting coupling, 2 x 2 liter
Mixing canulas, pink, 100 pcs.

Accessories	Item no.
Foot-operated switch	511136
Conversion kit for foot-operated switch	172630
Set material cask with quick fitting coupling, 2 x 2 liter	172640
Stand small for duplicating silicones, 2 x 6 kg-cask	172631
Stand tall for duplicating silicones, 2 x 25 kg-cask	172632
Connecting tubes cpl., 2 pcs.	172633
Adapter-set for duplicating silicones, 2 x 1 kg-cask	172634
Adapter-set for duplicating silicones, 2 x 6 kg-cask	172636
Mixing canulas 1:1, inner Ø 5.4 mm, pink, 100 pieces	622002
Consumable materials	
Duplicating silicone, 1:1	
WaSil soft, hardness Shore A: 16-18, 2 x 1 kg, component A / B	656060
WaSil soft, hardness Shore A: 16-18, 2 x 6 kg, component A / B	656061
WaSil soft, hardness Shore A: 16-18, 2 x 25 kg, component A / B	656062
WaSil speed, hardness Shore A: > 24, 2 x 1 kg, component A / B	656070
WaSil speed, hardness Shore A: > 24, 2 x 6 kg, component A / B	656071
WaSil speed, hardness Shore A: > 24, 2 x 25 kg, component A / B	656072
Sidomix	
Item no. 172984	
Voltage	220–240 V / 50/60 Hz
Power consumption	1.6 A
Output	260 W
W x H x D with bracket	290 x 310 x 450 mm
W x H x D without bracket	200 x 310 x 440 mm
Weight	11.5 kg
Sound level	≤ 70 dB (A)