


**Wamix-Touch:** The high-capacity and maintenance-free vacuum pump and the speed regulation of the mixer unit makes the unit well-suited to all kinds of plaster, embedding compounds and silicones. Depending on the application three vacuum levels are adjustable. Since it just requires a single power supply connection, the selection of the location is very flexible.


**Wamix-Touch Injector:** Vacuum is achieved rapidly with this process even for large volumes, thus saving time. Another advantage is the robustness of the highly efficient, oil-free and maintenance-free vacuum generator. Compressed air supply is required.



Wall fastening, mounting set included in delivery

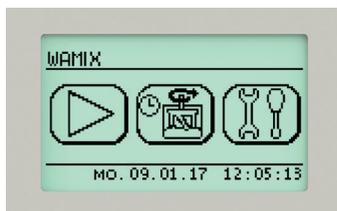

Touch screen for fast and easy operation

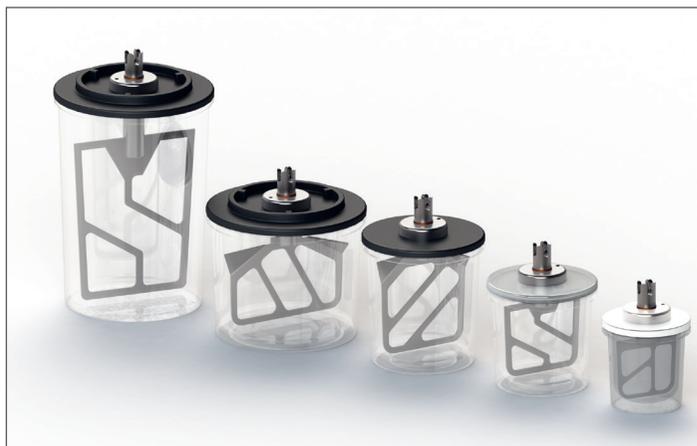

1000 ml, 600 ml, 350 ml, 200 ml and 75 ml

## MADE IN GERMANY

Wassermann is a medium-sized, owner-operated company with head office in Hamburg. We design and produce high-quality technical devices in Germany. Our main users are dental laboratories, dental clinics and various industrial areas worldwide.

### Why a Wassermann product?

Wassermann products are manufactured from high-quality materials and the proportion of plastics used is intentionally kept as low as possible. Our devices are known for their suitability for daily use, functionality and longevity. As a result, the users place high levels of trust in the Wassermann brand. We feel it is our duty to support our customers on a daily basis when it comes to safety, service and productivity. Technological perfection, optimum quality and especially sustainability are the benchmarks we set. Thanks to the combination of these values, the family-owned company, since its establishment in 1927, has gained worldwide recognition as a manufacturer of high-quality products.



Further information are available on our website at [www.wassermann.hamburg](http://www.wassermann.hamburg)



/WassermannDentalmaschinen

Your Dental Dealer

## Wamix-Touch and Wamix-Touch Injector



## Vacuum Mixers




Wassermann Dental-Maschinen GmbH  
Rudorffweg 15-17 · 21031 Hamburg · Germany  
Phone: +49 (0)40 730 926-20/ -24 · Fax: +49 (0)40 730 37 24  
export@wassermann-dental.com · www.wassermann.hamburg  
www.facebook.com/WassermannDentalmaschinen

# Vacuum Mixers

## Wamix-Touch & Wamix-Touch Injector

The **Wamix**-mixers by Wassermann impress with their perfect and bubble free results. They are easy to operate and extremely robust, which makes work much easier for sure. A touch screen and the newly developed control panel allow the user to set the relevant parameters quickly and easily using self-explanatory symbols. The technological highlights of the Wamix Series provide the foundation for first class products - quick and reproducible.

**Innovative:** The vacuum and speed values set by the user are operated by a microcontroller and are automatically readjusted.

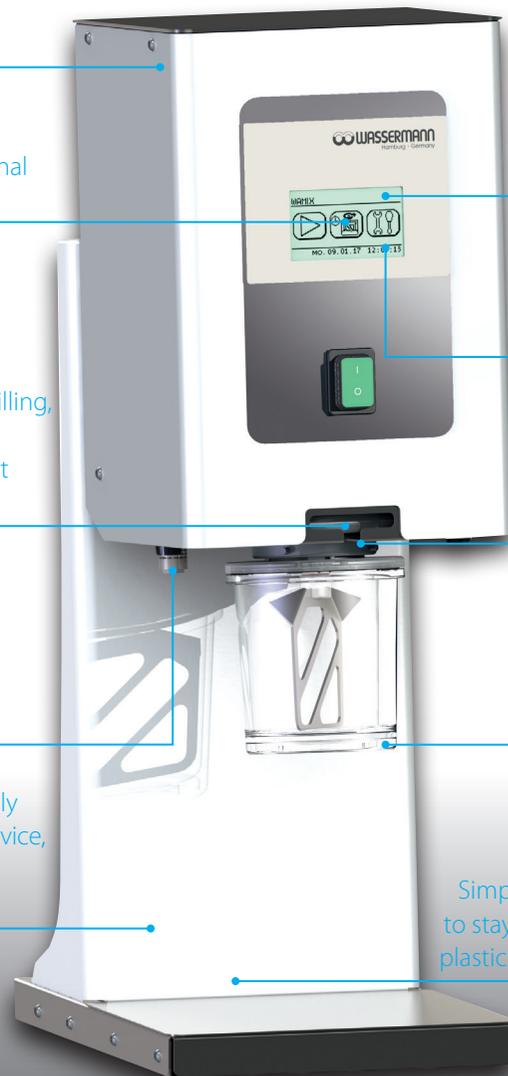

High quality metal components, almost maintenance-free

Visual and audible signal at programme end

Visual safety system helps to identify overfilling, additionally with protective filter against taking up liquids

Silencer Wamix-Touch Injector (Wamix-Touch: Silencer integrated in the unit)

Can be installed flexibly as a wall or counter device, table stand available as extras

Simple cleaning due to stay-clean surfaces, plastic powder coated

Quick and reproducible processes: You can save 30 programs with your choice of names and access them at any time

Beaker vented only after confirmation, post-vacuum also reduces gas inclusions

Inclined motor for easy insertion of the mixer beaker, comfortable single-hand operation, optimum mixing results

Automatic mixing when beaker is applied, programmable pre-mixing mode with right and left operation, adjustable change of rotation direction

### Major Features

- Reproducibility at the push of a button: Vacuum and speed monitoring guarantee homogeneous mixing of all kinds of plaster, embedding materials and silicones
- **Included in delivery:**  
Beaker, medium cpl. 350 ml  
Mounting set for wall fastening  
Compressed air hose, 2 m **Wamix-Touch Injector**

| Accessories   | Item no.                    |
|---|-----------------------------|
| Table stand (W 208 x H 450 x D 275 mm)                | 171744                      |
| Beaker, extra small cpl. 75 ml                        | 171295                      |
| Beaker, small cpl. 200 ml                             | 171230                      |
| Beaker, medium cpl. 350 ml                            | 171283                      |
| Beaker, large cpl. 600 ml                             | 171237                      |
| Beaker, extra large cpl. 1000 ml                      | 171269                      |
| Vacuum washing beaker, 600 ml                         | 171275                      |
| Compressed air hose per meter                         | 320026                      |
| <b>Wamix-Touch Injector</b>                           |                             |
| <b>Wamix-Touch</b>                                    | <b>Item no. 171970</b>      |
| <b>Wamix-Touch Injector</b>                           | <b>Item no. 171971</b>      |
| Voltage   | 100–240 V / 50/60 Hz        |
| Power consumption max. <b>Wamix-Touch</b>             | 0.44 A (Si 4.0 A)           |
| Power consumption max <b>Wamix-Touch Injector</b>     | 0.17 A (Si 4.0 A)           |
| Output <b>Wamix-Touch</b>                             | 94 W                        |
| Output <b>Wamix-Touch Injector</b>                    | 40 W                        |
| Speed (Motor control with intelligent regulation)     | 100 – 550 min <sup>-1</sup> |
| Vacuum <b>Wamix-Touch</b>                             | approx. 120 mbar (absolute) |
| Vacuum <b>Wamix-Touch Injector</b>                    | approx. 100 mbar (absolute) |
| Mixing time   | 999 sec. adjustable         |
| W x H* x D (wall fastening)                           | 200 x 310 x 200 mm          |
| W x H x D (with stand)                                | 210 x 590 x 300 mm          |
| Weight <b>Wamix-Touch</b>                             | 9.3 kg                      |
| Weight <b>Wamix-Touch Injector</b>                    | 7.5 kg                      |
| Sound level   | ≤ 80 dB (A)                 |
| Safety class  | IP 22                       |
| Compressed air connection <b>Wamix-Touch Injector</b> | min. 5 bar                  |

\*Leave at least 300-350 mm free space beneath the unit.

Other voltages on request · Technical changes reserved · 2021/06