

● MODEL TRIMMER HSS-88

The classic among wet trimmers



Technical information

Powerful, flexible and durable: There are good reasons why the versatile HSS-88 is Wassermann's most purchased trimmer. Its powerful motor and special starter mean 1000 Watts of nearly wear- and service-free operation at the highest cutting speeds. The HSS-88 is especially versatile and can use all Wassermann discs and orthodontic attachments. An outer housing cools the motor, so there is no risk of dirt entering. And of course it goes without saying that the trimmer operates with very little noise and no vibration.

- versatile, powerful wet trimmer
- can be used with all Wassermann discs
- high speed - 3000 min⁻¹
- very powerful, nearly maintenance-free motor, splash protected according to IP 54
- sturdy, oxidation-free bench with longitudinal slots and side openings – no water spills
- supplied with solenoid valve for automatic water supply
- powerful disc and internal rinsing
- easy to open using rear handle
- large work section 14.5 x 10.5 cm

● Technical data and accessories

HSS-88 without grinding disc	Item no.: 168997
Voltage*	220–240 V / 50/60 Hz
Power consumption	4.6 A
Output	1000 W
W x H x D	330 x 335 x 410 mm
Weight	19.8 kg
Sound level	≤ 70 dB (A)
Speed	3000 min ⁻¹
Duty	S6 - 40 %
Table dimensions projecting	120 mm

HSS-88 incl. grinding disc	Item no.:
incl. diamond grinding disc	168998
incl. Dia-Quick grinding disc (partly diamonded)	168996
incl. corundum grinding disc	168999
incl. foil grinding disc	168991


Incl. in delivery	Item no.:
Inclined support	168037
Supply hose 1 m	320006
Discharge hose 1 m	320017
Accessories	Item no.:
Conversion kit for HSS-88, Dia-Quick and foil grinding disc	168071
Grinding device KMS according to B. E. G. G.	168080
Directional light GTE	168081
Template for upper jaw to trim the lower jaw tube level	168082
Wall console	168056
Sealing grease for trimmer, 25 ml	629002
Truing apparatus to roughen corundum grinding disc	710001
Spare diamond grinding disc incl. puller and grease	168063
Spare Dia-Quick grinding disc (partly diamonded), incl. grease	162063
Spare corundum grinding disc, incl. puller and grease	168024
Spare foil grinding disc, 1 set = 4 pieces, incl. grease	162013
Spare foil grinding disc, fine, 1 set = 4 pieces, incl. grease	162023
Puller for grinding disc change	164026

*Other voltages on request

Technical changes reserved



HSS-88 with special accessories for orthodontics

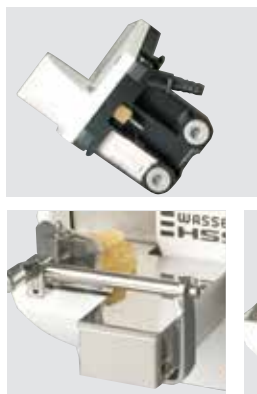

HSS-88 with inclined support

Major features at a glance

- most-purchased trimmer from Wassermann
- can be used with all Wassermann discs and orthodontic attachments
- 3000 min⁻¹
- powerful 1000 Watt motor
- multiple expansion options
- the grinding attachment KMS can be delivered also on swiss- and EOS-version

www.wassermann-dental.com