

● **DOSING AND MIXING UNIT**
SIDOMIX

A head for the right mixture


Technical information

Perfect mixture at the press of a button: with the Sidomix, mixing and dosing of low-viscous, addition cross-linking duplicating silicones is now child's' play. Because, with the simple press of a button the desired components are mixed 1:1 and dosed correctly- fully automatically, totally homogenous and fully bubble free. The innovative mixing head is the secret: its spatially separated outlet openings prevent the components from mixing - even when attaching and removing the mixing cannulas. Practical in its application: the containers for the A and B components can be attached very simply to the unit with the quick-couplings.

- automatic mixing and dosing unit for low-viscous, addition cross-linking duplicating silicones
- components are mixed and dosed 1:1 at the push of a button, homogenously and bubble free
- material and time saving with loss free filling of several cells with one mixing cannula
- practical adapter with quick coupling for simple attachment of 1 kg and 6 kg casks, support frame available for 6 and 25 kg casks
- foot switch available as option
- fillable material containers with quick fitting coupling for components A and B included in delivery scope

● Technical data and accessories

Sidomix	Item no.: 172984
Voltage*	220–240 V / 50/60 Hz
Power consumption	1.6 A
Output	260 W
W x H x D incl. bracket	290 x 310 x 450 mm
W x H x D without bracket	200 x 310 x 440 mm
Weight	11.5 kg
Sound level	≤ 70 dB (A)



Consumable materials Duplicating silicone, 1:1

	Item no.:
WaSil soft*, hardness Shore A: 16-18, WaSil speed**, hardness Shore A: > 24, 2 x 1 kg, component A / B	656060* / 656070**
WaSil soft*, hardness Shore A: 16-18, WaSil speed**, hardness Shore A: > 24, 2 x 6 kg, component A / B	656061* / 656071**
WaSil soft*, hardness Shore A: 16-18, WaSil speed**, hardness Shore A: > 24, 2 x 25 kg, component A / B	656062* / 656072**

*Other voltages on request

Incl. in delivery	Item no.:
Set material cask incl. quick fitting coupling, 2 x 2 liter	172640
Mixing canulas, pink, 100 pcs.	622002
Accessories	Item no.:
Foot-operated switch	174175
Conversion kit for foot-operated switch	172630
Set material cask with quick fitting coupling, 2 x 2 liter	172640
Stand small for duplicating silicones, 2 x 6 kg-cask	172631
Stand tall for duplicating silicones, 2 x 25 kg-cask	172632
Connecting tubes cpl., 2 pcs.	172633
Adapter-set for duplicating silicones, 2 x 1 kg-cask	172634
Adapter-set for duplicating silicones, 2 x 6 kg-cask	172636
Mixing canulas, pink, 100 pcs.	622002

Technical changes reserved

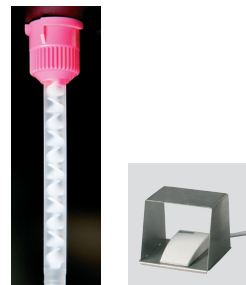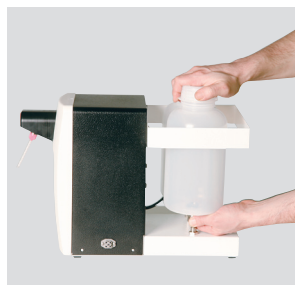

Major features at a glance

- specially recommended for use of WaSil-duplicating silicones
- automatic mixing and dosing unit for duplicating silicones
- absolutely homogeneous and bubble-free mixture
- mixing head with separate outlet openings – makes mixing of components impossible
- practical adapter with quick fitting coupling for silicone casks A and B
- high-quality manufacture, housing plastic-powder-coated
- foot-operated switch available as option

www.wassermann-dental.com