


## • HYDRAULIC PRESS

### WW-33

The correct pressure day after day

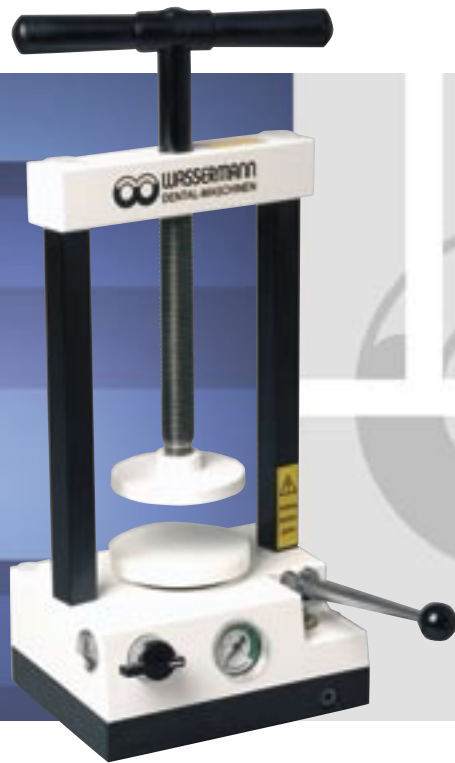

### Technical information

Safe, infinitely adjustable and incredibly strong: Finally, the correct pressure for your flasks with the WW 33 hydraulic press. Designed for daily use, with a working pressure that can be adjusted up to 6 tonnes – for the correct pressing of 1–3 or 2–4 flasks. And so efficiently, thanks to the valve which can be opened to allow the piston to return to its starting position automatically – ready again for operation without losing time. Of course, the oxidation-free piston is made of super-strong material, as is the indestructible design in drawn steel. The Wassermann credo: Safety from tradition!


- break-resistant hydraulic press in drawn steel for pressing 1–3 or 2–4 flasks
- pressure infinitely adjustable up to 6 tonnes
- oxidation-free piston in super-strong material
- automatic piston return after valve opening
- quality pressure gauge with precision display recessed in housing
- support separation of 152 mm to allow use of larger flasks and flask clamps
- high stability, integrated oil level display

## ● Technical data and accessories

|                          |  |
|--------------------------|--|
| <b>WW-33, 2–4 flasks</b> | <b>Item no.: 126998</b>                          |
| <b>WW-33, 1–3 flasks</b> | <b>Item no.: 126999*</b>                         |
| WxD/ D incl. pump lever  | 200 x 158 / 295 mm                               |
| Distance between columns | 152 mm   |
| Distance plunger plate   | max. 272 / min. 102 mm<br>max. 212 / min. 42 mm* |
| Weight                   | 16.0 kg, 14.7 kg*                                |
| Working pressure         | up to 6 tons                                     |
| Working table            | Ø 95 mm  |

| Accessories                                | Item no.: |
|--|-----------|
| Universal flask                            | 129990    |
| Special steel tension bracket for 3 flasks | 129992    |
| Special steel tension bracket for 2 flasks | 129993    |
| Special steel tension bracket for 1 flask  | 129994    |
| Universal flask bracket for 3 flasks       | 129995    |
| Universal flask bracket for 2 flasks       | 129996    |
| Universal flask bracket for 1 flask        | 129997    |
| Hydraulic oil 500 ccm                      | 629005    |

Technical changes reserved

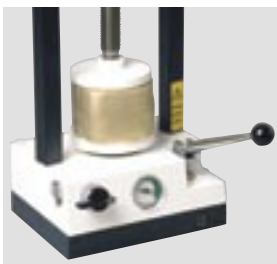

WW-33 1-3 flasks




All-steel press 1-3 flasks

### All-steel press 1–3 flasks

**Item no.: 129999**

|                          |                               |
|--------------------------|-------------------------------|
| W x max. H / H x D       | 210 x max. 565 / 395 x 100 mm |
| Distance between columns | 155 mm                        |
| Distance plunger plate   | max. 218 / min. 47 mm         |
| Weight                   | 9.5 kg                        |
| Working pressure         | up to 1 ton                   |

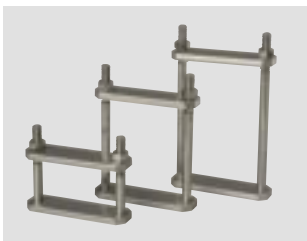

Special steel tension bracket



Universal flask bracket



## Major features at a glance WW-33

- fracture-resistant hydraulic press is made of drawn steel with oxidation-free piston
- working pressure infinitely adjustable up to 6 tonnes
- automatic piston return
- suitable for larger flasks and flask clamps
- housing plastic-powder-coated

[www.wassermann-dental.com](http://www.wassermann-dental.com)