

● WAX SCALDING UNIT

WAPO-Ex 12 II

Effective and efficient
boiling out



Technical information

The Wapo-Ex 12 combines competence, many years of experience and the latest technology to produce an outstanding product. This well-designed free standing unit with integrated spray gallows scalds up to 12 flask halves simultaneously and fully automatically. The innovative multifunction touch-screen terminal offers diverse functions and reproducible procedures by simple programming. The Wapo-Ex 12 is very durable and stands out thanks to its flexible usage, programmable timer and automatic wax removal. The electronic water level gauge, the leakage signal and the inlet limitation ensure additional safety.


- free standing unit for scalding up to 12 flask halves
- multifunction touch-screen terminal for easiest operation
- efficient usage thanks to its programmable timer and automatic wax removal
- stainless steel construction
- safety features with visual and audible signals
- easy to clean

● Technical data and accessories

Wapo-Ex 12 II	Item no.: 170970	Incl. in delivery	Item no.:
Voltage*	220-240 V / 50/60 Hz	Spray gallows	830145
Power consumption	12 A	Service set (lime remover and wax remover)	170330
Output	max. 3300 W	Supply hose with water inlet valve and gasket	830011
W x H (without spray gallows) x D	650 x 900 x 620 mm	Discharge hose with hose band clamp Ø 12-20 x 9	320061/322006
W x H (with spray gallows) x D	650 x 1600 x 620 mm	Spray arm tube with seal	830160
Weight	87 kg	2 flask baskets	2 x 830205
Sound level	≤ 70 dB (A)	2 covers	830230 / 830235
Water capacity	approx. 50 l	Accessories	Item no.:
Water feed rate	4.5 l / min	Service set (lime remover and wax remover)	170330
Water drain height	470 mm		

*Other voltages on request

Technical changes reserved

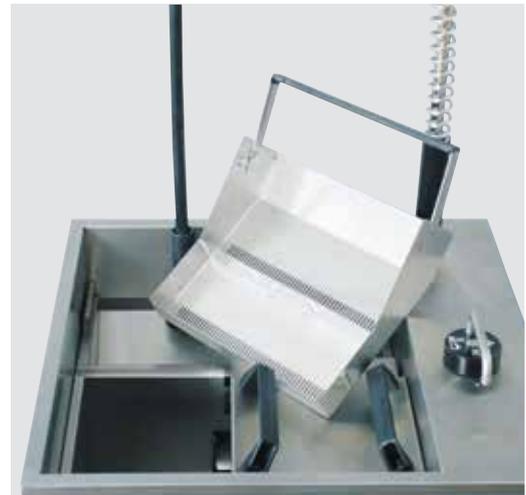

Major features at a glance

- scalding of up to 12 flask halves
- easy to operate
- very durable and energy saving because of the intelligent heating control and the high quality insulation
- stainless steel housing
- efficient usage thanks to its programmable timer and automatic wax removal

www.wassermann-dental.com