

● **WAX SCALDING UNIT****WAPO-Ex 8 III**

Compact shape, huge performance

**Technical information**

Compact, mature, durable: The Wapo-Ex 8 is the ideal solution if space is limited but the demand is great. Up to 8 flask halves are scalded simultaneously and fully automatically. The sophisticated multifunction touch-screen terminal allows diverse functions and reproducible procedures to be created through the simplest programming. First-class workmanship, very high-quality and maintenance-free components, plus an electronic water level gauge, guarantee maximum safety, low servicing costs and a long service life. User-friendly: The spray protector and integrated spray head for manual scalding on the recessed filter cover.

- compact tabletop unit for scalding up to 8 flask halves
- informative multifunction touch-screen terminal for easy operation
- thermostat-protected, gasketless immersion pump in bronze; maintenance-free
- energy-saving tank thermal insulation
- spray system set into side for easy access to flask baskets
- all stainless steel construction
- visual and audible water-level indicator
- easy to clean without tools

## ● Technical data and accessories


Wapo-Ex 8 III	Item no.: 170981
Voltage*	220-240 V / 50/60 Hz
Power consumption	9 A
Output	2000 W
W x H x D	610 x 540 x 420 mm
Weight	42.0 kg
Sound level	≤ 70 dB (A)
Water capacity	18.5 l

Incl. in delivery	Item no.:
Service set (filters, lime remover, wax remover)	170021
PVC hose 19 x 4 (3/4") – 1.5 m	320030
Hose clamp Ø 20 - 32 x 12.2	322017

Accessories	Item no.:
Service set (filters, lime remover, wax remover)	170021
Mobile stand	170120
PVC hose 19 x 4 (3/4") – per meter	320030

\*Other voltages on request

Technical changes reserved

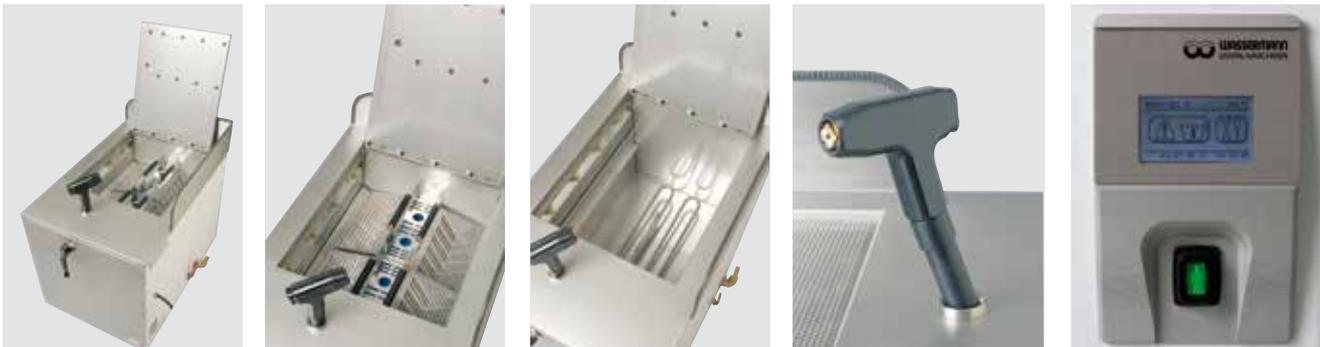

### Major features at a glance

- compact scalding of up to 8 flask halves
- multifunction touch-screen terminal
- easy to operate
- durable and energy-saving
- stainless steel housing
- integrated hand spray and spray protector

[www.wassermann-dental.com](http://www.wassermann-dental.com)