


# Trimmer incl. Grinding Discs

Device alternatives	Item no.
HSS-99 incl. diamond grinding disc	164998
HSS-88 incl. diamond grinding disc	168998
HSS-88 incl. Dia-Quick grinding disc (partially diamond coated)	168996
HSS-88 incl. corundum grinding disc	168999
HSS-88 incl. foil grinding disc	168991
HSS-AZ incl. Dia-Quick grinding disc (partially diamond coated)	162997
HSS-AZ incl. foil grinding disc	162996
HSS-ZA incl. Dia-Quick grinding disc (partially diamond coated)	161996
HSS-ZA incl. foil grinding disc	161997



HSS-99, with diamond grinding disc and orthodontic accessories

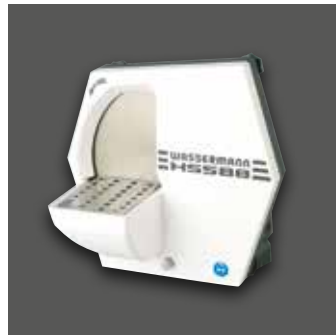

HSS-88, can be equipped with all Wassermann discs



HSS-AZ, designed for perfect grinding in the dental clinics



HSS-ZA, in the rear with polishing spindle

## MADE IN GERMANY


Wassermann is a medium-sized, owner-operated company with head office in Hamburg. We design and produce high-quality technical devices in Germany. Our main users are dental laboratories, dental clinics and various industrial areas worldwide.

### Why a Wassermann product?

Wassermann products are manufactured from high-quality materials and the proportion of plastics used is intentionally kept as low as possible. Our devices are known for their suitability for daily use, functionality and longevity. As a result, the users place high levels of trust in the Wassermann brand. We feel it is our duty to support our customers on a daily basis when it comes to safety, service and productivity. Technological perfection, optimum quality and especially sustainability are the benchmarks we set. Thanks to the combination of these values, the family-owned company, since its establishment in 1927, has gained worldwide recognition as a manufacturer of high-quality products.



## Grinding Discs

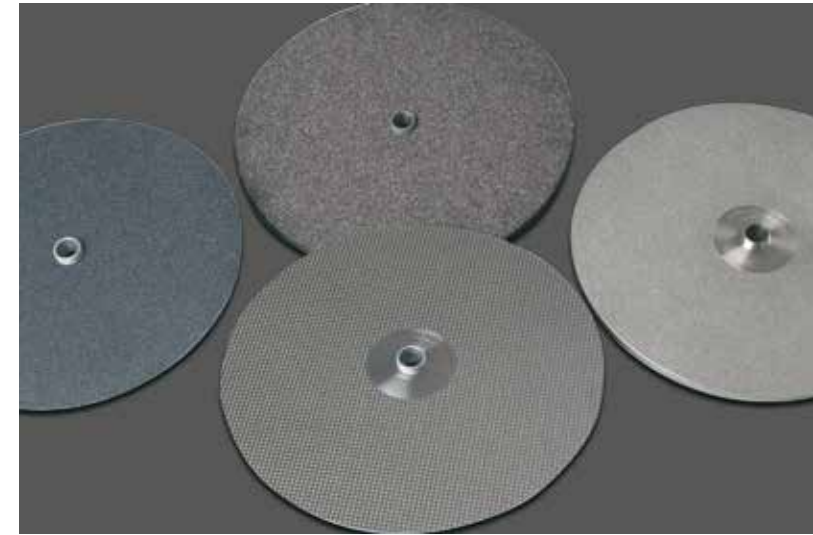

For HSS-99, HSS-88, HSS-AZ, HSS-ZA



Further information are available on our website at [www.wassermann.hamburg](http://www.wassermann.hamburg)

Your Dental Dealer



Wassermann Dental-Maschinen GmbH  
 Rudorffweg 15-17 · 21031 Hamburg · Germany  
 Phone: +49 (0)40 730 926-0 · Fax: +49 (0)40 730 37 24  
[info@wassermann-dental.com](mailto:info@wassermann-dental.com) · [www.wassermann.hamburg](http://www.wassermann.hamburg)


# Grinding Discs

Wassermann discs set standards. Available for every plaster hardness and offering the user top quality to meet the highest demands. From the versatile foil grinding disc to Wassermann's own diamond grinding discs – our experience is your benefit.



## Diamond Grinding Disc for HSS-99 and HSS-88

- Wassermann diamond coated discs make easy work of even the hardest plaster – minimal grinding pressure, maximum time saving
- Fully diamond coated for smooth surfaces, model surfaces are always even despite the aggressive grinding power
- Rough and fine grinding in one process
- Extremely long service life – diamonds support each other over the whole surface, hence they cannot break off
- Excellent use for orthodontics

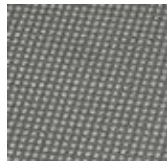

## Dia-Quick Grinding Disc for HSS-88 and HSS-AZ/-ZA

- Partially diamond coated, island shape design
- Rapid results even with low grinding pressure
- Long service life
- Maximum safety due to electronically balanced supporting disc



## Corundum Grinding Disc for HSS-88

- Approx. 10 mm thick silicon carbide coating
- Self-regenerating therefore constant and powerful grinding over a long period
- Complete surface glued to supporting disc for maximum safety
- Electronically balanced, therefore no vibration



## Foil Grinding Disc for HSS-88 and HSS-AZ/-ZA

- Silicon carbide coating – usable on both sides
- Flexible, stable design
- Electronically balanced supporting disc
- Disc not glued to supporting disc, hence easily to be replaced
- Fast and easy replacement of worn discs

## Major Features

### Included in delivery:

Sets of spare grinding discs	Item no.
Spare diamond grinding disc for <b>HSS-99</b> , fully diamond coated, incl. sealing grease and puller for grinding disc change	164063
Spare diamond grinding disc for <b>HSS-88</b> , fully diamond coated, incl. sealing grease and puller for grinding disc change	168063
Spare Dia-Quick grinding disc for <b>HSS-88, HSS-AZ, HSS-ZA</b> , partially diamond coated, incl. sealing grease	162063
Spare corundum grinding disc for <b>HSS-88</b> , 10 mm thick, incl. sealing grease and puller for grinding disc change	168024
Spare foil grinding disc for <b>HSS-88, HSS-AZ, HSS-ZA</b> , 1 set = 4 pieces, incl. sealing grease	162013
Spare foil grinding disc, fine, for <b>HSS-88, HSS-AZ, HSS-ZA</b> , 1 set = 4 pieces, incl. sealing grease	162023
Conversion kit for <b>HSS-88</b> for Dia-Quick- and foil grinding disc, incl. sealing grease and puller for grinding disc change	168071