

MADE IN GERMANY

Wassermann is a medium-sized, owner-operated company with head office in Hamburg. We design and produce high-quality technical devices in Germany. Our main users are dental laboratories, dental clinics and various industrial areas worldwide.

Why a Wassermann product?


Wassermann products are manufactured from high-quality materials and the proportion of plastics used is intentionally kept as low as possible. Our devices are known for their suitability for daily use, functionality and longevity. As a result, the users place high levels of trust in the Wassermann brand. We feel it is our duty to support our customers on a daily basis when it comes to safety, service and productivity. Technological perfection, optimum quality and especially sustainability are the benchmarks we set. Thanks to the combination of these values, the family-owned company, since its establishment in 1927, has gained worldwide recognition as a manufacturer of high-quality products.



WP-Ex 10 II



Polishing Compact Unit

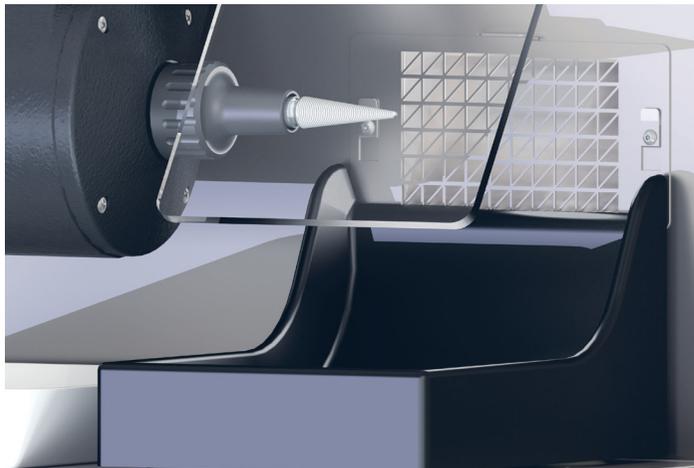

WP-Ex 10 II with special accessory suction protective grid



WP-Ex 10 II with special accessory daylight LED spot



Further information are available on our website at www.wassermann.hamburg



/WassermannDentalmaschinen


Your Dental Dealer _____



Wassermann Dental-Maschinen GmbH
 Rudorffweg 15-17 · 21031 Hamburg · Germany
 Phone: +49 (0)40 730 926-20/ -24 · Fax: +49 (0)40 730 37 24
 export@wassermann-dental.com · www.wassermann.hamburg
 www.facebook.com/WassermannDentalmaschinen

Polishing Compact Unit WP-Ex 10 II

The WP-Ex 10 II is Wassermann's most compact polishing unit for first class polishing and suctioning in the smallest space - a high quality construction, made almost entirely of stainless steel. In addition to its powerful, longlife and maintenance-free squirrel-cage motor it also convinces with a powerful and low noise suction.

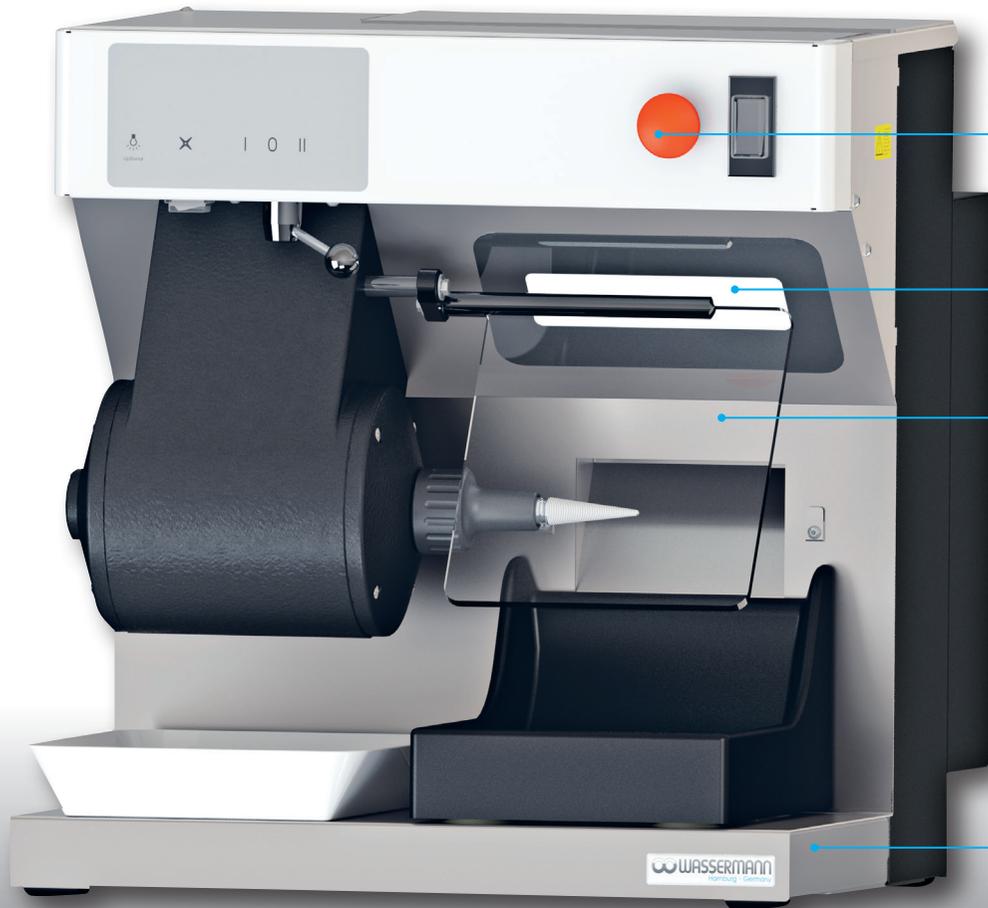

Emergency off-switch and restart protection guarantee safety at work

Integrated daylight LEDs

Impact-resistant, adjustable safety screen

Work area made of stainless steel

Major Features

- Compact and powerful polishing unit for daily use at laboratories even in the smallest space
- Motor with jacket-cooling, dust protected
- 2 settings, 1500 resp. 3000 min⁻¹
- Micro fine dust filterbag with large capacity which is easy to replace
- Low-noise extraction
- Very convenient: The easy to clean soft-rubber polishing trough can be removed very quickly
- Daylight LED spot available as special accessory
- **Included in delivery:**
Brush spindle, 92 mm
Micro fine dust filterbag
Storage bowl (W 257 x H 40 x D 198 mm)
Polishing trough
Supply line

Accessories	Item no.
Daylight LED spot	146390
Suction protective grid, 1 piece	146396
Brush spindle right, 92 mm	152022
Brush spindle right, 125 mm	152043
Stone holder right	152056
Drill chuck right, 1-6 mm clamping	152049
Storage bowl (W 257 x H 40 x D 198 mm)	830409
Spare micro fine dust filterbag, 1 set = 2 pieces	611005
WP-Ex 10 II	Item no. 146971
Voltage	220–240 V / 50/60 Hz
Power consumption	2.2 A
Output	500 W
Fan capacity	445 m ³ /h (free-blowing)
W x H x D	455 x 485 x 505 mm
Weight with motor	37.0 kg
Sound level	≤ 65 dB (A)
Speed	1500 / 3000 min ⁻¹